manual

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Item: 40772

EN

NO

SV


www.rubicson.com
Box 50435 Malmö
Sweden
2021-09-02


Rubicson®

Overview

Massage gun

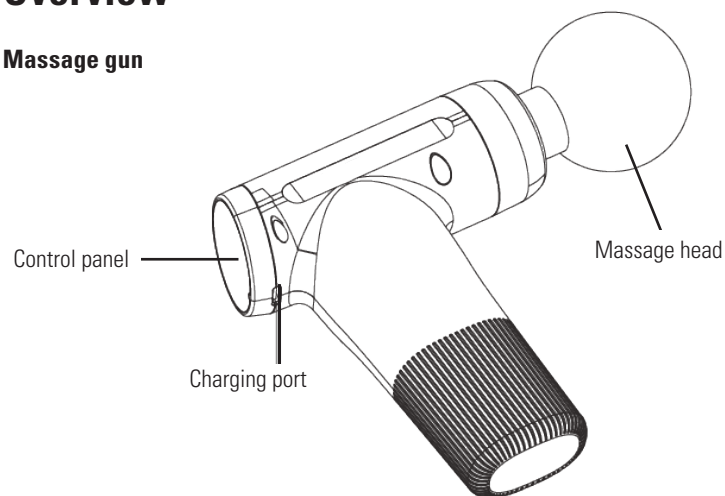
Massage heads

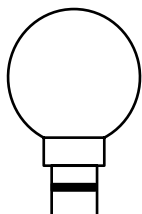
Spherical head

For soft tissue and sensitive muscle groups.

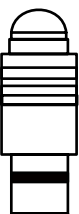
Cone head

For deep tissue. Treats trigger points and massages joints.

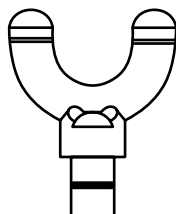U-shaped head

Used for e.g. neck and heel tendons. Use around the muscles rather than straight on.

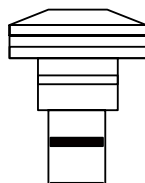
Flat head

For large muscle groups and flat areas such as the back, chest and legs.

Control panel

1. Gear indicators
2. Power button
3. Charging indicator

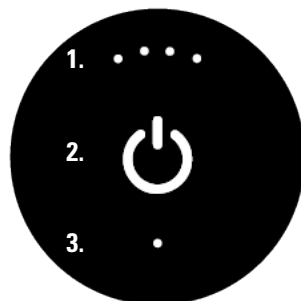Specifications

Speeds:	4, up to 3200 rpm
Battery:	1x rechargeable 8.4 V battery (2x 2000 mAh)
USB-C input:	5 V, 2 A
Charging time:	3 hrs
Run time:	3-4.5 hrs
Motor power:	30 W
Noise level:	60-73 dB

Use

1. Insert a massage head into the massage gun.
2. Press and hold the power button on the control panel for two seconds to turn the massage gun on.
3. Press the power button again to change the speed. The gear indicators on the control panel show the current speed.
4. Place the massage gun on the area you want to massage and gently press and move it around. Note! Do not use the massage gun for more than 15 consecutive minutes and do not massage the same body part for more than 4-5 minutes. Do not massage directly on your skin, make sure there are clothes or a dry towel in between.
5. Press and hold the power button for two seconds to turn the massage gun off. Remove the massage head by pulling it out.

Charging

The battery indicator flashes red when the battery level is low.

Connect the USB-C cable to the charging port of the massage gun and to a USB adapter. The power indicator keeps flashing red during charging. The massage gun is fully charged when the indicator turns green.

Safety information

Do not massage directly on your skin. Make sure there is a garment or a dry towel between you and the massage gun.

This appliance is not intended for use by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they are supervised or have received instructions concerning the use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.

High speed and pressure must only be used on soft tissue, and not for the head or bones.

Frequent use of high speed and massage with high pressure on the same part may cause bruises. Stop using the massage gun if it feels uncomfortable or painful.

Keep the massage gun away from high heat or low temperatures.

Do not disassemble or repair the product yourself.

The massage gun is for adults only.

Do not leave the massage gun unsupervised during charging or if it is powered on.


Disconnect the charging cable from the massage gun when it is fully charged. Excessive charging may cause damage to the battery.

Keep the massage gun away from water. It is not waterproof.

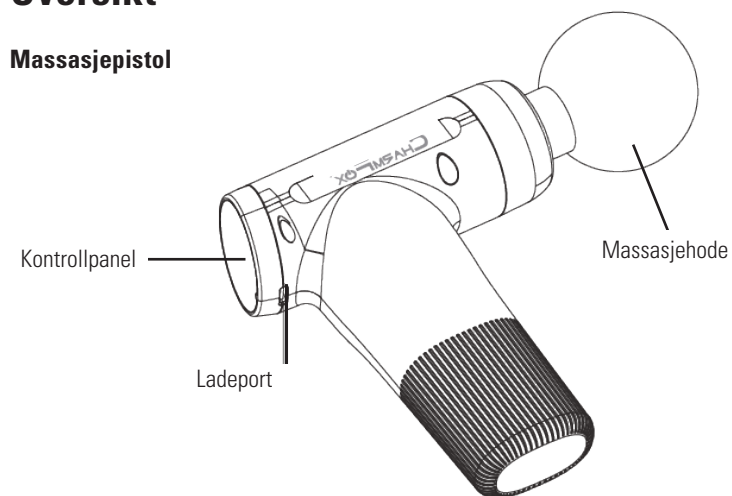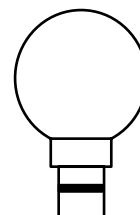
Do not use the massage gun if it looks damaged or works abnormally.

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