manual

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Art no: 52390

EN

NO

SV


www.rubicson.com
Box 50435,
Malmö, Sweden
2024-06-20


Rubicson®

Overview

1. Massage head piston
2. Charging port
3. Speed level/battery indicator
4. Decrease speed level
5. Power on/off
6. Increase speed level

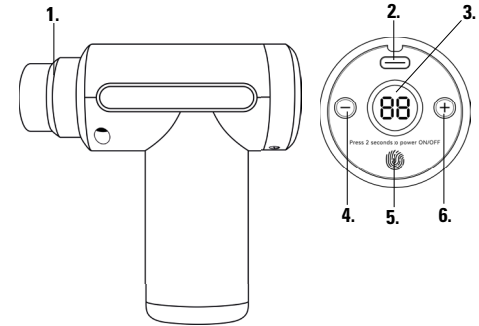Massage heads

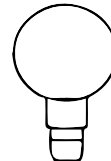Ball head

Large muscle groups, back, thighs, chest.

Forked head

Suitable for the neck and spine

Bullet head

Deep muscle tissue, joints, soles.

Flat head

Universal use.

Specifications

Speeds:	1–20, 1800–2800 rpm
Battery:	Rechargeable, 7.4 V, 2000 mAh
Power input:	5 V/1 A (charger sold separately)
Charging time:	Approx. 4 h
Run time:	3.5–4 h
Size and weight:	130x133x50 mm, 560 g (excl. massage head)
Included:	Massage gun, USB charging cable (45 cm), 4 massage heads, manual

Before first use

Charge the massage gun fully before using it for the first time.

Use

Charging

Connect the included USB charging cable between the charging port **(2)** and a USB power adapter (sold separately). The ring around the battery indicator **(3)** starts flashing. The indicator shows “99” when the massage gun is fully charged.

Note! Do not use the massage gun during charging.

Mounting a massage head

Attach the massage head you want to use to the massage gun piston **(1)** by pressing the bottom of the head into the front of the massage gun. Make sure the head is securely fastened before moving to the next step.

Power on/off and speed settings

Press and hold your finger for 2 seconds on the area marked Power on/off touch control **(5)**. The massage gun battery level indicator **(3)** lights up and shows the remaining battery capacity (from 99 % to 20 %). Press the head of the massage gun to the muscles you want to massage.

Press the speed level buttons **(4, 6)** to increase and decrease the speed level, the display shows the speed level. It is recommended to start with a low setting and slowly increase the speed for a more enjoyable massage. There are 20 speed levels available where level 1 runs at 1800 rpm and level 20 runs at 2800 rpm.

Press and hold your finger for 2 seconds on the Power on touch control **(5)** to turn off the massage gun.

Note! The massage gun automatically stops after 15 minutes of continuous use or if too much force is used during a massage. Allow the massage gun to rest for at least 1 minute before using it again. Press the button to increase speed level **(6)** to restart the massage. To ensure the longevity of the motor, turn the machine off after 1 hour of use and let it rest for 30 minutes before using the massage gun again.

Maintenance

Cleaning

Clean the massage gun and the massage heads using a damp towel and dry them with a soft cloth. Do not use any detergents on any of the parts.

Safety information

This product contains a rechargeable Li-ion battery that cannot be replaced.

Do not store the product with the battery completely discharged for extended periods of time. To avoid damaging the battery, charge the product once every 3 months during extended periods of storage.

This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.

Children shall not play with the appliance. User maintenance shall not be made by children without supervision.

Do not use the massage gun on injured body parts before consulting a doctor.

Never attempt to disassemble the product. Repairs should only be performed by a qualified technician. By repairing it yourself there is a risk of electric shock and/or short circuit.

Keep the device away from moisture, water, high temperatures, and direct and strong sunlight.

Do not use the massage gun to directly massage dry skin to avoid injuries. The massage area should either be separated from the massage gun by a dry towel or clothing or used with massage oil. It is not recommended to massage a single muscle or part of the body for more than 60 seconds at a time.

Only use high speeds and high-pressure settings on soft tissue parts of the body to avoid causing pain and discomfort. Do not use high speed and high pressure on the head and any bone parts of the body. Frequent use of high speeds and high pressure on the same body part may cause abrasions. If you feel pain and/or discomfort, stop using the massage gun immediately.


Keep fingers, hair and/or any other part of the body away from the moving telescopic rod of the massage gun and the heat dissipation vent of the motor to prevent injuries.

Do not cover the heat dissipation vent or insert any foreign objects into the motor.

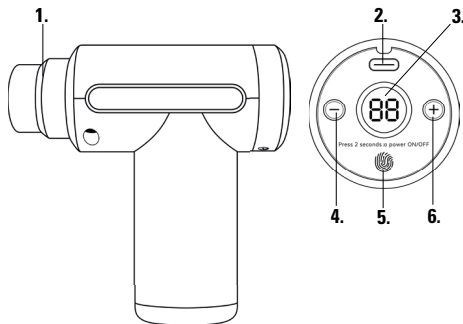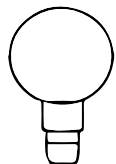
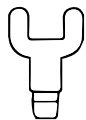
Do not leave the massage gun running unsupervised.

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